



Universal Couchtop™

SBRT Enabled | DoseMatch™ Technology



CIVCO
MEDICAL SOLUTIONS




DoseMatch™ Technology

Assists in dose accuracy between planning and treatment

Easy - Identical indexing makes it easy to transfer patient setups from simulation to treatment.

Efficient - Using a matching overlay reduces time during treatment planning and assists in providing more accurate treatments.

Compatible - CT overlays are compatible with a variety of models from GE, Philips, Siemens and Toshiba.

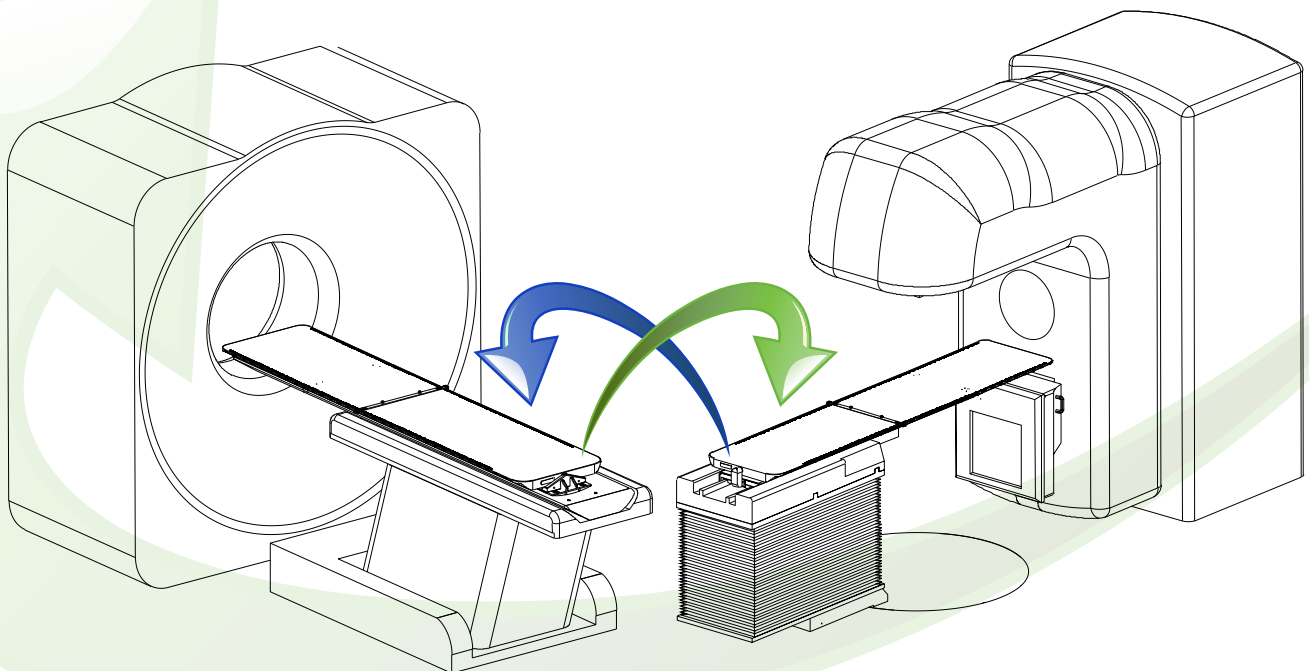

Check with your account representative to confirm DoseMatch Technology is available for your CT Simulator.

SBRT Enabled

Comfortable immobilization for hyperfractionation

Unique - The Wing Profile of the Universal Couchtop allows for easy attachment of Body Pro-Lok™ bridges via a small adapter, providing anterior immobilization and breath restriction options.

Custom - Customize your setup to each individual patient. Body Pro-Lok Bridges easily adjust to accommodate a variety of patient sizes and utilize your existing two-pin indexable products.

Variety - Body Pro-Lok accessory options include bridges, forehead restraint, shoulder restraint, Clam-Lok™ cushion, respiratory plate and respiratory belt.

IGRT Friendly

Low attenuation and lightweight

Metal Free - Metal-free extension seam reduces artifacts while maintaining a strong attachment for extensions.

Visibility - The Universal Couchtop's homogeneous design provides high quality images for IGRT.

Accuracy - Compatible with CIVCO's Protura™ 6DOF Robotic System for increased setup accuracy.



Additional Features

- Self-Locking extensions
- Utilizes Prodigy™ 2 indexing at 7cm increments providing more setup flexibility
- Gray inferior surface for better visualization of SSD and field light
- Optional accessory rail and patient handles
- SBRT and Respiratory Restriction packages for convenient accessory ordering



Available Accessories



Type-S™ Extension




Rectangular Extension



Posifix® Extension (Available soon)



Varian End Plate

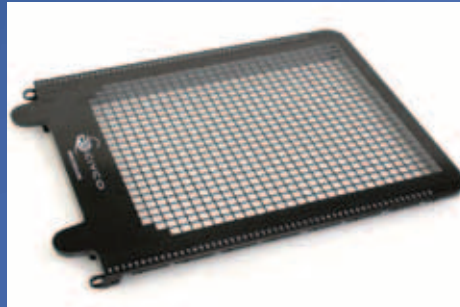

Palliative Extension - Grid



Palliative Extension - Solid



Prodigy™ 2 Lok-Bar™



Accessory Rail



Patient Handle



Universal Couchtop™

SBRT Enabled | DoseMatch™ Technology

800.842.8688 | +1 712.737.8688 | WWW.CIVCO.COM

COPYRIGHT © 2009, CIVCO AND POSIFIX ARE REGISTERED TRADEMARKS OF CIVCO MEDICAL SOLUTIONS. UNIVERSAL COUCHTOP, DOSEMATCH, TYPE-S, PRODIGY, LOK-BAR, BODY PRO-LOK, CLAIM-LOK AND PROTURA ARE TRADEMARKS OF CIVCO. ALL OTHER TRADEMARKS ARE PROPERTY OF THEIR RESPECTIVE OWNER. ALL PRODUCTS MAY NOT BE LICENSED IN ACCORDANCE WITH CANADIAN LAW. 2009P0527 Rev.A

